

# CLAF



## Introduction

The EscO CLAF™ Ceiling Laminar Airflow is a containment equipment module utilized to provide a zoned laminar down-flow airflow. The purpose of this equipment is to provide laminar airflow to a working area for product protection.

Ceiling Laminar Airflows are customizable units that are:

- Utilized as open restricted access barriers over filling and capping machines.
- Stand alone units mounted via eye bolts and drop rods over specific applications.
- Stand alone units mounted over mobile stands for mobile aseptic zones.

## Basic Principles

- Room air is drawn pre-filtered via an EU6 prefilter before entering via the perforated diffuser into supply plenum.
- A special baffle system that channels the airflow through the gel-sealed HEPA filters as downflow supply creating an aseptic work zone with low noise.

## Key Features

- Easy to clean.
- HEPA/ULPA gel-sealed design better than conventional gasket sealed.
- Sentinel silver microprocessor control with audio/visual alarms for downflow velocity.
- Zoned magnehelic gauges for filter loading.
- Energy efficient teardrop lightings away from downflow.
- Emergency stop

## Optional Features

- Remote mounted main control panel
- Splashproof electrical outlets
- PVC Curtains
- Splashproof electrical outlets

## Guide to Models

### CLAF-PC-1      -AC-H13

Product Code	Materials of Construction	Electrical Code	External Width	External Dimension	External Height	Fan	Filter Type
CLAF	PC - Powder coated EG Steel	1 - 220-240 VAC 50/60 Hz	TBD*	TBD*	TBD*	AC	H13
	S1 - Stainless Steel 304	2 - 110-130 VAC 50/60 Hz					H14
	S2 - Stainless Steel 316 L	3 - 100-110 VAC 50/60 Hz					U15
							Others

\* To be discussed according to customers' preference


## Electrical Supply Options

1 = 220-240 VAC 50 / 60 Hz

2 = 110-130 VAC 50 / 60 Hz

3 = 100-110 VAC 50 / 60 Hz

# ESCO GLOBAL NETWORK

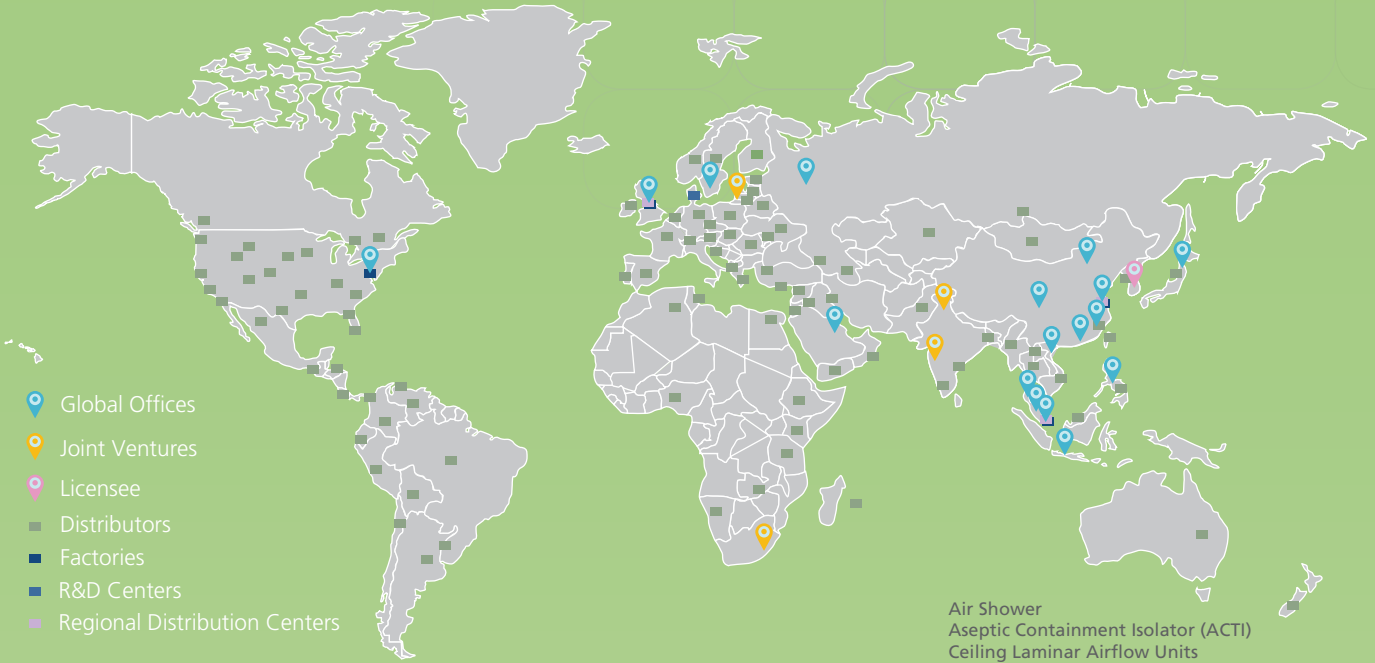

- Air Shower
- Aseptic Containment Isolator (ACTI)
- Ceiling Laminar Airflow Units
- Cleanroom Transfer Hatch
- Containment Barrier Isolator (CBI)
- Compounding Aseptic Isolator
- Compounding Aseptic Containment Isolator
- Downflow Booth (DFB)
- Dynamic Floor Label Hatch
- Dynamic Pass Box
- Evidence Drying Cabinet
- Laminar Flow Storage Cabinet
- General Processing Platform Isolator (GPPI)
- Healthcare Platform isolator
- Laminar Flow Horizontal Trolley
- Laminar Flow Straddle Units, Single and Double Laminar Flow Vertical Trolley
- Pass Box
- Soft Wall Cleanroom
- Sputum Booth
- Ventilated Balance Enclosure (VBE)
- Weighing and Dispensing Containment Isolator (WDCI)

Since 1978, Esco has emerged as a leader in the development of controlled environment, laboratory and pharmaceutical equipment solutions. Products sold in more than 100 countries include biological safety cabinets, fume hoods, ductless fume hoods, laminar flow clean benches, animal containment workstations, cytotoxic cabinets, hospital pharmacy isolators, and PCR cabinets and instrumentation. With the most extensive product line in the industry, Esco has passed more tests, in more languages, for more certifications, throughout more countries than any biosafety cabinet manufacturer in the world. Esco remains dedicated to delivering innovative solutions for the clinical, life science, research and industrial laboratory community. [www.escoglobal.com](http://www.escoglobal.com).



**Esco Pharma Pte Ltd**  
 21 Changi South Street 1 Singapore 486777 • Tel: +65 65420833  
 Fax: +65 65426920 • Email: [csis.pharma@escoglobal.com](mailto:csis.pharma@escoglobal.com)

**Esco Technologies, Inc.**  
 2512 Metropolitan Drive, Suite 120 B  
 Feasterville- Trevose, PA 19053-6738  
 Tel: 215 322 2155 • Email: [eti.pharma@escoglobal.com](mailto:eti.pharma@escoglobal.com)

**Esco Gb Ltd**  
 Unit 2 R-Evolution @ Gateway 36, Kestrel Way, Barnsley, S70 5SZ  
 Tel: +44 (0) 1226 360 799 • Email: [egb.pharma@escoglobal.com](mailto:egb.pharma@escoglobal.com)

**Esco Global Offices:** Bahrain | Bangladesh | China | India | Indonesia | Italy | Japan  
 Malaysia | Philippines | Russia | Singapore | South Africa | South Korea | Thailand  
 United Kingdom | USA | Vietnam



CLAE\_sellsheet\_AA\_v8\_12112017  
 Esco can accept no responsibility for possible errors in catalogues, brochures and other printed materials.  
 Esco reserves the right to alter its products and specifications without notice. All trademarks and logos in this  
 material are the property of Esco and the respective companies.